

## Specifications:

**Hearit**  
HÖRAT

Communicator

Hearit MT



**description** Hearit MT is an effective and uncomplicated communicator with a built in telecoil and an acknowledged superior sound quality.

**applications** Hearit MT is intended for people with reduced hearing who want a very easy-to-handle, robust and powerful solution.

**function** Hearit MT is very easy to use. It has one large control for on/off and volume adjustments and an illuminated switch for microphone or telecoil selection. It's built in telecoil is for use in rooms with a teleloop system. With all controls at maximum the total amplification is 56 dB (50 dB at 1,6 KHz). The AGC safeguards the ears against damage from sudden sounds that are too loud such as being dropped on a table. Hearit MT takes readily available standard AA batteries, which are easy to install.


### Technical specifications

User controls:	On/Off switch Volume control M/T switch
Preset control:	Bass cut
Dimensions:	150 x 40 x 30 mm
Weight with batteries:	110 g
Battery type:	Two 1,5V penlight, LR6, AA
Battery life:	About 220 hours
Power consumption:	6-9 mA
Max. gain (acoustic):	56 dB
Max. output level (acoustic):	120 dB SPL
Impedance of headphones/neckloop:	16-32 Ω

### Accessories




890-88-629-00  
standard headphones

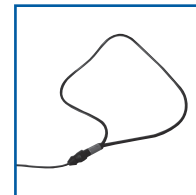

890-88-604-00  
neckloop



890-15-141-00  
conference  
microphone

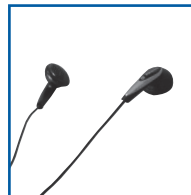

890-88-601-00  
ear buds



890-88-616-00  
headphones on the



895-88-044-00  
bag for Hearit MT and  
headset

### Standard package 890-88-530-00



Hearit MT



standard  
headphones



AA batteries  
2 pcs are



bag for Hearit  
MT and  
headset

### Multiple package for TV 890-88-620-00



Connecting cable  
to TV box  
2 x 1,5 m



890-88-622-00  
TV box to  
connect up to four  
Hearit to the TV



392-88-107-00  
scart cable to  
connect TV box  
to TV



392-88-103-00  
scart cable to  
connect Hearit  
to TV (8 m)

### TV Cable

**Phonic Ear®**

Kongebakken 9 • 2765 Smørum

Tel.: **+45 39177101** • Fax: **+45 39277900** • E-mail: **mail@phonicear.com** • Web: **www.phonicear.com**

Hearit MT Spec sheet Rev. B 0911